

Enteral Feeding Pump



KL-5031N



KL-5051N
with automatic flush function



KL-5021A

KellyMed | Jev&Kev
SINCE 1994



One-stop Solution for Enteral Nutrition

Beijing KellyMed Co., Ltd.

ADD: 6R International Metro Center, No. 3 Shilipu, Chaoyang District, Beijing, 100025, China

TEL: +86-182 1056 5691 FAX: +86-10-6558 7908

WEB: www.kelly-med.com EMAIL: africa@kelly-med.com


Enteral Feeding Set for Single Use.

1.Features

- 1)Our dual-layer co-extrusion tubes use TOTM (DEHP free) as plasticizer. The inner layer doesn't contain colorant. The purple color of the outer layer can prevent misusing with IV sets.
- 2)Compatible with various feeding pumps and liquid nutrition containers.
- 3)Its internationally ENFit® connector can be used for various nasogastric feeding tubes. Its ENFit® connector design can prevent feeding tubes from accidentally fitting into IV sets.
- 4)Its ENFit® connector is very convenient for feeding medicine and flushing tubes.
- 5)We have different models and specifications to meet different clinic needs.

2.Types


Our products can be used for nasogastric feeding tubes, nasogastric stomach tubes, enteral nutrition catheter and feeding pumps.

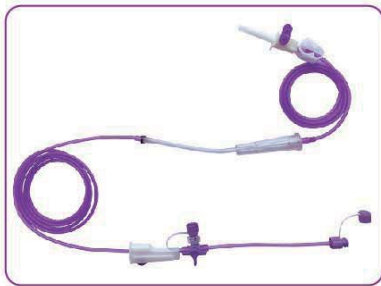


Spike Set



Bag Set

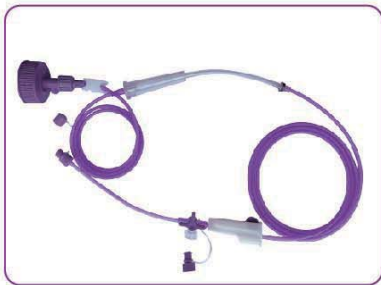

Screw Cap Set

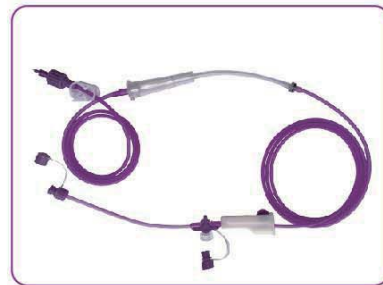

Screw Spike Set

The standard length of silicone tube is 11cm and 21cm.
11 cm is suitable for roller feeding pump, and 21cm is suitable for peristaltic feeding pump.
We can customize any length for the feeding tubes according to customers' requirement.

3.Models and specifications

The bag set has 500ml, 800ml, 1000ml, 1200ml, 1500ml, 2000ml.

		Model	Package
Pump Feeding		Spike Set	JP2-1-101, JP2-1-102, JP2-1-103, JP2-1-104, JP2-1-105, JP2-1-106
		Bag Set	JP2-2-101, JP2-2-102, JP2-2-103, JP2-2-104, JP2-2-105, JP2-2-106
	Screw Set	Screw Cap Set	JP2-3-101, JP2-3-102, JP2-3-103, JP2-3-104, JP2-3-105, JP2-3-106
		Screw Spike Set	JP2-3-107, JP2-3-108, JP2-3-109, JP2-3-110, JP2-3-111, JP2-3-112
Gravity Feeding		Spike Set	JP2-1-001, JP2-1-002
		Bag Set	JP2-2-001, JP2-2-002
	Screw Set	Screw Cap Set	JP2-3-001, JP2-3-002
		Screw Spike Set	JP2-3-003, JP2-3-004

4. The connection for Enteral Feeding Sets and Nutrition Container

Model	Picture	Nutrition Container	How to connect
Spike Set	With air inlet	With rubber bottle stopper	Insert the spike to the bottle stopper
Set for Bottle	Three-in-one Connector	Plastic wide-mouth bottle	Special three-in-one design enables to connect with bottle mouth directly by screwing it
		Glass bottle	Special three-in-one design enables to connect with bottle mouth directly
Set for bag	Screw connector	Bag nutrition container	Screw the threaded connector to connect with the bag directly.